

# VersaVision

## Digital Specimen Radiography System

### Introducing VersaVision from Faxitron®

Featuring the highest resolution detector in the industry, VersaVision is ideal for imaging core biopsies and a broad range of specimens. VersaVision provides instant high-quality images for immediate sample verification, and its compact size makes it ideal for benchtop placement in radiology or wherever you need easy access to imaging.



### VersaVision FEATURES

#### Highest Resolution Imaging

With a 21 lp/mm detector that provides up to 5X geometric magnification, VersaVision is ideal for imaging core biopsies and a broad range of specimens at over 60 lp/mm.

#### Immediate Verification of Core Samples

High-quality images are instantaneously displayed on a high-resolution 24" monitor, allowing for instant sample verification, helping to reduce procedure time and improve workflow.

#### Compact and Versatile

The lightweight, fully shielded x-ray cabinet can be placed directly on the bench top in radiology or wherever you need easy access to imaging.

#### Easy to Use

With no additional training or specialized requirements needed for x-ray shielding, VersaVision plugs into any A/C outlet and is powered on with a single button.

#### Advanced Automatic Exposure

With the click of a mouse, the system automatically selects the appropriate exposure time and kV settings for the sample.

#### Ideal for a Range of Samples

The 5X mag tray provides extremely high-resolution (60+ lp/mm) for the smallest samples, and the laser centering guide ensures proper sample positioning.

#### Reliable and Secure Data Transfer

The VersaVision system can send multiple annotated images to PACS at a touch of a button.


#### Single Click Calibration

Auto object detection allows for simple, error-free calibration.

#### Easy to Maintain

The stainless steel inner chamber has been designed for easy cleaning and maintenance.

## All Top 10 US Cancer Centers Have Multiple Faxitron Systems

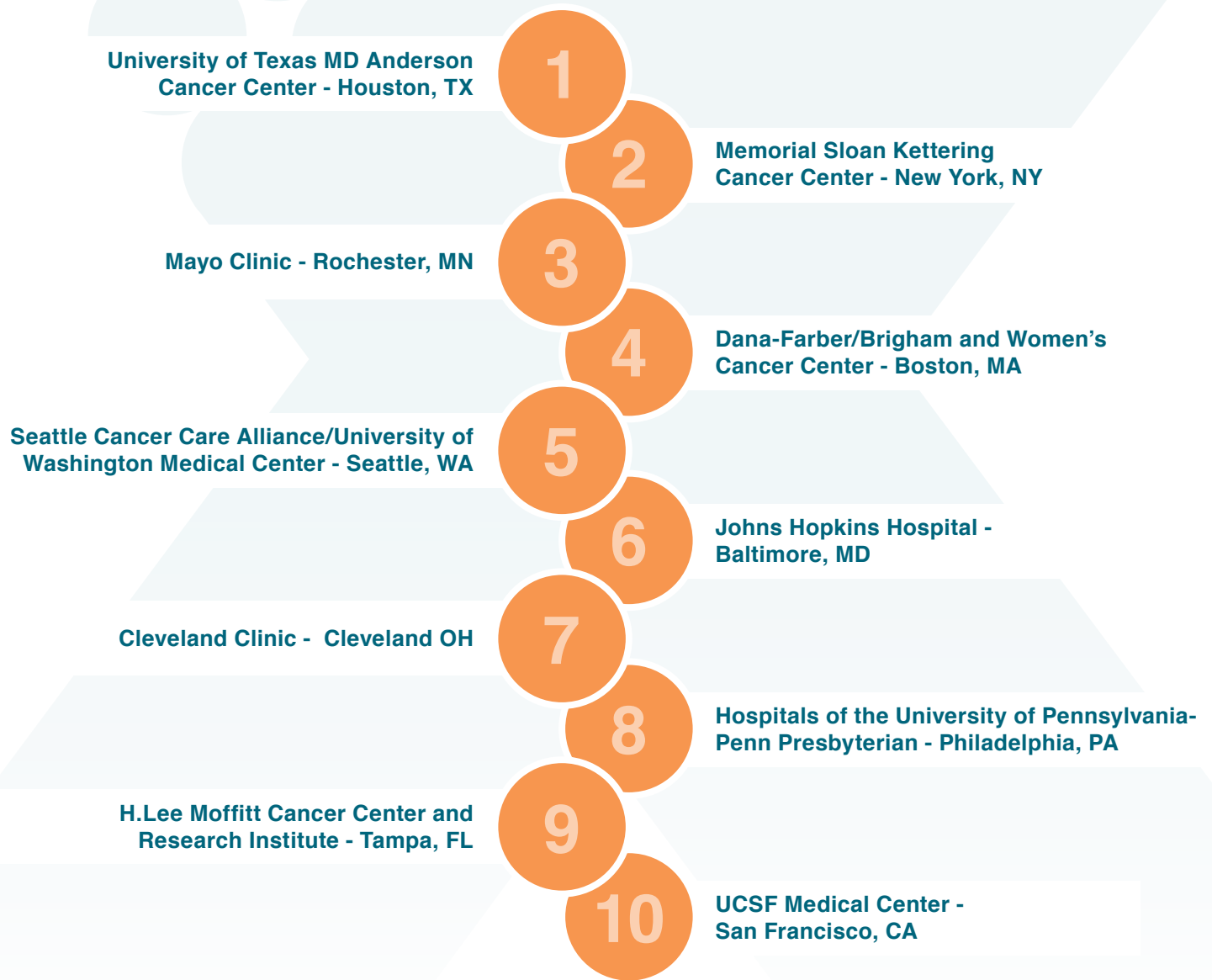

Source Data from US News and World Report - Top Cancer Centers 2017-18. <http://health.usnews.com/best-hospitals/rankings/cancer>

### FAXITRON BIOPTICS, LLC

3440 E Britannia Drive, Suite 150,  
Tucson, Arizona USA 85706

General inquiries: +1.520.399.8180

Toll Free: +1.877.910.0030

Fax: +1.520.399.8182

[sales@faxitron.com](mailto:sales@faxitron.com)

[www.faxitron.com](http://www.faxitron.com)